The Perils of Multiplex Gastrointestinal Pathogen Panels: Pseudo-outbreaks of Salmonellae and *Entamoeba histolytica* in Immunocompromised Hosts

Elizabeth Robilotti, MD, MPH;^{1,2,4} Eleanor Powell, PhD;³ Shauna Usiak, MPH;² Ying Taur, MD MPH;^{1,4} N. Esther Babady, PhD;^{1,3} Mini Kamboj, MD^{1,2,4}

Two distinct clusters of gastroenteritis due to Salmonellae and *Entamoeba histolytica* (EH) were identified using a multiplex gastrointestinal pathogen panel (GPP) at a tertiary-care cancer center. Despite temporo-spatial overlap, our investigation did not corroborate transmission or true infection. In clinical practice, GPPs may render false-positive results.

Infect Control Hosp Epidemiol 2018;39:867-870

Acute diarrheal illness, common among hospitalized immunocompromised patients has many etiologies, including infection, medications (eg, laxatives), and chemotherapy side effects.¹ Rapid assessment of infection is essential to ensure timely therapeutic intervention and implementation of isolation precautions in vulnerable patient populations. Multiplex PCR-based panels for gastrointestinal pathogens offer rapid and simultaneous analysis of multiple bacterial, viral, and parasitic pathogens, with highly convenient same-day turnaround of results. Additionally, these assays are less dependent on the technical expertise of laboratory personnel than on traditional methods. For the 2 predominant multitarget gastrointestinal pathogen panels (GPPs) cleared by the US Food and Drug Administration (FDA), sensitivity and specificity are high. However, clinical experience with routine implementation, particularly in specialized populations and for pathogens uncommon in high-income countries, is insufficient.²⁻⁴ This report describes the investigation of 2 clusters of uncommon infections, one Entamoeba histolytica (EH) and the other Salmonellae, in hospitalized patients. Upon in-depth investigation, both clusters were determined to be pseudo-outbreaks related to GPP testing.

METHODS

The study was conducted at the Memorial Sloan Kettering Cancer Center, a 471-bed tertiary-care cancer hospital in New York, New York, where Luminex xTAG GPP (Luminex Diagnostics, Toronto, Canada) had been implemented in 2014 for the evaluation of stool from patients with diarrheal illness. The infection control department performs routine surveillance for all pathogens included on the panel. The Luminex GPP has been cleared by the FDA for 14 targets: 8 bacteria, 3 viruses, and 3 parasites. Among them, 5 are reportable to the New York State Department of Health, including Salmonellae and EH. Additionally, under New York State (NYS) and New York City (NYC) Health Code laws, laboratories are required to submit certain isolates for confirmatory testing. According to the protocol, isolate recovery for all samples positive by GPP for Salmonellae were attempted using Hektoen Enteric agar plates (Hardy Diagnostics, Santa Maria, CA). All positive raw stool samples on GPP are routinely saved and frozen at -80°C.

Cluster 1

From April 19, 2017, to May 4, 2017, 5 cases of EH were identified via routine laboratory evaluation of clinical cases of non-travel-related diarrhea using the Luminex xTAG GPP. Among them, 3 were adult hematopoietic stem cell transplant (HSCT) recipients; 2 cases were diagnosed 2 days apart on the same HSCT unit, while the third occurred in an HSCT recipient residing on another unit. The 2 additional cases occurred children with recent inpatient admission or exposure to our facility's pediatric day hospital.

Cluster 2

Between July 1 and July 5, 2017, 3 cases of infection with *Salmonella* spp were identified by GPP in hospitalized patients admitted to the same oncology unit. The first 2 cases occurred within 1 day in patients residing in adjoining rooms. The third case was identified 3 days later on the same unit. The 3 patients were cared for by different clinical teams with very limited overlap of support staff.

No clinical or ancillary staff reported contemporaneous diarrheal illness at the time of either cluster, and no common diet- or procedure-related links could be established among the cases. All patients received targeted antibacterial or antiparasitic therapy for the identified pathogen. Clinical and demographic characteristics of all cases are described below (Table 1).

RESULTS

Testing

The mean overall inpatient and outpatient GPP testing volume per month for 2017 was 317, compared to 311 for 2016 (Figure 1). This difference was not significant (P=.60; 95% CI, -33.8 to 21). No selective increase in testing was detected among the units with cases. Notably, the 3 cases of *Salmonella* identified on GPP in October 2017 were identified in outpatients without known prior temporo-spatial overlap. One of

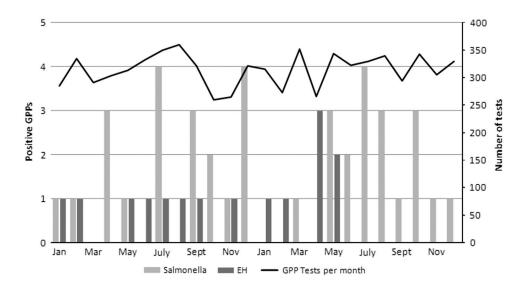


FIGURE 1. Salmonellae- and Entamoeba histolytica-positive GPPs and overall tests per month January 2016 through December 2017.

these patients remained GPP positive in December. No copathogens were identified on GPP or by other routine clinical diagnostic testing (Table 1).

Retesting

A second analysis of the positive specimens was performed with another FDA-approved multiplex assay, FilmArray (BioFire Diagnostics, Salt Lake City, UT). Residual samples, when available, were also sent to the New York City and/or New York State Department of Health for additional confirmatory testing, which included EH and *E. dispar* polymerase chain reaction (PCR; suspected EH cases) and stool culture (suspected Salmonellae cases). These alternative methods did not confirm EH or *Salmonella* spp in any of the respective specimens. Confirmatory testing results are summarized in Table 1.

DISCUSSION

The clusters described here were initially compelling because possible outbreaks because of corroborating gastrointestinal symptoms and close temporo-spatial proximity among cases. Clinical symptoms of both *Salmonella* and EH are nonspecific and can mimic other types of colitis common among oncology patients such as *Clostridium difficile* colitis, cytomegalovirus colitis, and gastrointestinal graft-versus-host disease. Further investigation revealed no epidemiologic exposure for the identified pathogen among all patients. While confirmatory results were in progress, clinicians elected to treat patients due to their immunocompromised state and lack of alternate microbiologic explanation before ultimately attributing symptoms to medication side effects. While no adverse events related to treatment of suspected EH or *Salmonella* were identified, these patients did endure unnecessary exposure to antibacterials and amebicides. Additionally, substantial personnel time and resources at the local and public health levels were required for the investigation of these pseudooutbreaks. The associated costs may be an unanticipated expense for health systems that transition to multiplex testing from traditional methodologies.

Entamoeba histolytica, the causative agent of amebiasis, is a common parasitic infection transmitted to humans through contaminated food and water. While it is ubiquitous in many parts of the world, it remains an uncommon cause of diarrheal illness among hospitalized patients in developed countries. The clinical manifestations of intestinal EH infection can range from asymptomatic colonization to severe colitis. For the xTAG GPP, poor specificity for EH was recognized in an early evaluation of the assay for multipathogen testing.⁵

Salmonellae, gram-negative bacteria, frequently cause foodborne gastrointestinal illness. Disease manifestations range from self-limited diarrheal illness to bacteremia and invasive infections especially among immunocompromised hosts. Hospital-based outbreaks due to Salmonellae have been reported, and the occurrence of any cluster should be promptly investigated.^{6,7}

Multiplex panels have several distinct advantages, especially in terms of turnaround time from sample collection to actionable results and the limited need for specialized laboratory training. However, these advantages may be undercut by real-world performance in certain patient populations. Currently, data are limited in the reported literature on the real-world performance characteristics of GPP and their associated clinical implications in the evaluation of diarrheal illness in developed nations or among immunocompromised hosts. The specificity issue with Salmonellae and EH detection

		Underlying						Acid				Е.		
_	Age/	Disease	Signs/		_	Luminex	BioFire GPP		Giardia	/		dispar		EH
Case	Gender	(+t-HSCT)	Symptoms	Chemotherapy	Treatment	Result ^a	Result	Stain	FA	FA	Stain	PCR		PCR
Enta	moeba hi	istolytica												
1	42/M	AML (Y)	Diarrhea	Y	Metronidazole, parmomycin	+	-	-	-	-	-	-	-	Cluster 1
2	66/M	AML (Y)	Diarrhea, colitis on imaging	Y	Metronidazole, parmomycin	+	-	-	-	-	-	-	-	
3	65/M	DLBCL (Y)	Diarrhea	Y	Metronidazole, parmomycin	+	-	-	-	-	-	-	-	
4	8/M	ALL (N)	Diarrhea, abdominal pain, fever	Y	Metronidazole, parmomycin	+	-	ND	ND	ND	ND	ND	ND	
5	5/F	ALL (Y)	Diarrhea, abdominal pain	Y	Metronidazole, parmomycin	+	-	ND	ND	ND	ND	ND	ND	
			-			Initial Luminex GPPª	Luminex GPP (repeat) ^a	BioFire GPP	MSK culture ^b		NYC & NYS DOH culture			
Saln	ıonella						-							
1	25/M	GCT (Y)	Diarrhea, abdominal pain, fever	Y	β-lactam	+	-	-	1	No growth	No g	rowth		Cluster 2
2	68/F	Lung (N)	Diarrhea	Y	β-lactam	+	-	-	1	No growth	No g	rowth		
3	41/F	ATLL (Y)	No symptoms	Ν	β-lactam	+	-	+	1	No growth N		JD		

TABLE 1. Patient Demographics and Gastrointestinal Pathogen Confirmatory Testing

NOTE. HSCT, hematopeotic stem cell transplant; GPP, gastrointestinal pathogen panel; FA, fluorescent antibody; PCR, polymerase chain reaction; EH, *Entamoeba histolytica*; AML, acute myeloid leukemia; DLBCL, diffuse large B-cell lymphoma; GCT, germ cell tumor; ATLL, adult T-cell leukemia/lymphoma; MSK, Memorial Sloan Kettering; NYC & NYS DOH, New York City and New York State departments of health; ND, not done.

^aAll Luminex median fluorescence intensity results exceeded the manufacturer's thresholds for a positive result.

^bOverall Salmonellae recovery in culture from January 2016 to December 2017 was 27% (10 of 37 isolates).

was highlighted in a small report from Netherlands, in which 11 of 43 additional positives using xTAG GPP could not be confirmed, including 5 EH and 4 *Salmonella* cases.⁸ These results, coupled with those of our study, support the need for confirmatory testing in epidemiologically implausible circumstances. A large health-system analysis of GPP implementation noted the large number of positive results of uncertain clinical significance as the main detraction to widescale implementation of these testing strategies, while others have identified GPP use in hospitalized patients as a potential diagnostic stewardship target.^{9,10}

In summary, clusters of a single pathogen identified by multiplex GPP diagnostic assay require investigation to confirm the likelihood of a true outbreak, especially in the absence of convincing epidemiologic data for exposure to the identified pathogen. Further evaluation of GPP diagnostic assays is necessary to document the true sensitivity and specificity of these multiplex assays in real-world clinical settings, especially among specialized patient populations.

ACKNOWLEDGMENTS

Financial support: The Memorial Sloan Kettering Cancer Center supported this study (Core grant no. P30 CA008748).

Potential conflicts of interest: Dr. Babady acknowledges grant support from Luminex. All authors report no conflicts of interest relevant to this article.

Affiliations: 1. Infectious Diseases, Department of Medicine, Memorial Sloan Kettering Cancer Center, New York, New York; 2. Infection Control, Department of Medicine, Memorial Sloan Kettering Cancer Center, New York, New York; 3. Clinical Microbiology Service, Department of Laboratory Medicine, Memorial Sloan Kettering Cancer Center, New York, New York; 4. Department of Medicine, Weill Cornell Medical College, New York, New York.

Address correspondence to Elizabeth Robilotti, MD, MPH, 1275 York Avenue Box 9, New York, NY 10065 (robilote@mskcc.org).

Received January 24, 2018; accepted March 28, 2018; electronically published May 10, 2018

© 2018 by The Society for Healthcare Epidemiology of America. All rights reserved. 0899-823X/2018/3907-0016. DOI: 10.1017/ice.2018.95

REFERENCES

- Cox GJ, Matsui SM, Lo RS, et al. Etiology and outcome of diarrhea after marrow transplantation: a prospective study. *Gastroenterology* 1994;107:1398–1407.
- Halligan E, Edgeworth J, Bisnauthsing K, et al. Multiplex molecular testing for management of infectious gastroenteritis in a hospital setting: a comparative diagnostic and clinical utility study. *Clin Microbiol Infect* 2014;20:O460–O467.
- Khare R, Espy MJ, Cebelinski E, et al. Comparative evaluation of two commercial multiplex panels for detection of gastrointestinal pathogens by use of clinical stool specimens. *J Clin Microbiol* 2014;52:3667–3673.
- 4. Laude A, Valot S, Desoubeaux G, et al. Is real-time PCR-based diagnosis similar in performance to routine parasitological examination for the identification of *Giardia intestinalis*, *Cryptosporidium parvum/Cryptosporidium hominis* and *Entamoeba histolytica* from stool samples? Evaluation of a new commercial multiplex PCR assay and literature review. *Clin Microbiol Infect* 2016;22:e191–e198.
- Navidad JF, Griswold DJ, Gradus MS, Bhattacharyya S. Evaluation of Luminex xTAG gastrointestinal pathogen analyte-specific reagents for high-throughput, simultaneous detection of bacteria, viruses, and parasites of clinical and public health importance. *J Clin Microbiol* 2013;51:3018–3024.
- Smith AM, Mthanti MA, Haumann C, et al. Nosocomial outbreak of *Salmonella enterica* serovar Typhimurium primarily affecting a pediatric ward in South Africa in 2012. *J Clin Microbiol* 2014;52:627–631.
- 7. Tan SH, Lawler J, Foster K, et al. Nosocomial transmission of *Salmonella typhimurium* in renal transplant recipients. *J Hosp Infect* 2010;75:241–242.
- Wessels E, Rusman LG, van Bussel MJ, Claas EC. Added value of multiplex Luminex Gastrointestinal Pathogen Panel (xTAG(R) GPP) testing in the diagnosis of infectious gastroenteritis. *Clin Microbiol Infect* 2014;20:O182–O187.
- Freeman K, Mistry H, Tsertsvadze A, et al. Multiplex tests to identify gastrointestinal bacteria, viruses and parasites in people with suspected infectious gastroenteritis: a systematic review and economic analysis. *Health Technol Assess* 2017;21: 1–188.
- Hitchcock MM, Gomez CA, Banaei N. Low yield of FilmArray GI panel in hospitalized patients with diarrhea: an opportunity for diagnostic stewardship intervention. *J Clin Microbiol* 2018;56 pii e01558–17; doi: 10.1128/JCM.01558-17.